



Aria™ Low-Profile Brackets

A Comfortable Bracket Producing Optimum Results




Optimizing Precision, Comfort, and Control of the MIM-Manufactured *Aria* Bracket

Aria Low-Profile Brackets harmonize the enduring precepts of classic orthodontic design with state-of-the-art engineering. The result is a complete orthodontic appliance system that achieves the optimum clinical outcome in metal aesthetics, as well as the comfort that patients seek. The patented Horizontal Scribe Line (HSL) provides ideal visual cues for precise placement. The advanced Metal Injection Molding (MIM) technology maintains the precision and control of the design, allowing for the outcomes that uncompromising practitioners require. A finely-tuned appliance design for predictable, efficient treatment, and discerning finishes.

ANTI-ROTATION FOR THE ROTH RX IS BUILT INTO THE ARCHWIRE SLOT


By placing the anti-rotation into the slot, the mesial and distal sides of the bracket are presented perpendicular to the tooth. This eliminates the parallax that is created as a result of placing the anti-rotation in the base of the brackets.

The anti-rotation in the slot allows for the clinician to more accurately place the bracket

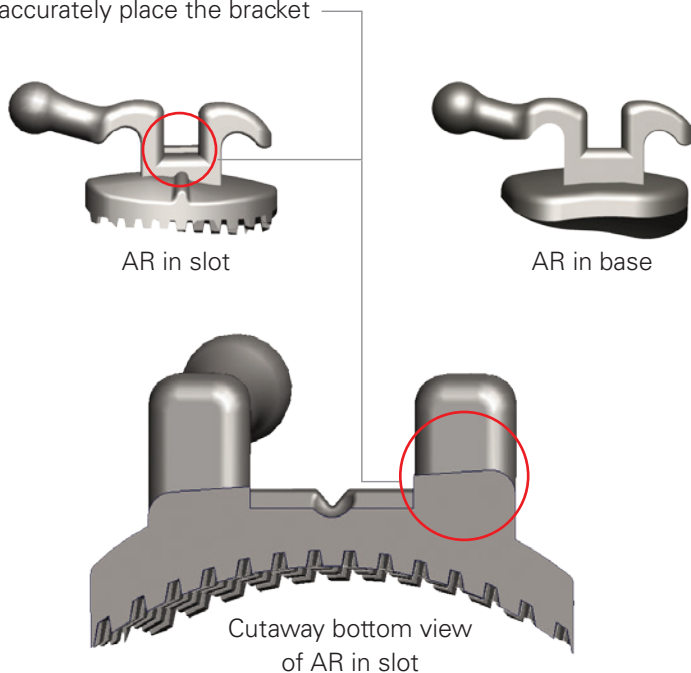

Rhomboid shape brackets and pads allow for maximum placement control

Patented Horizontal Scribe Line across the brackets and pads, mesial distal, creating additional visual clue, to enhance bracket placement

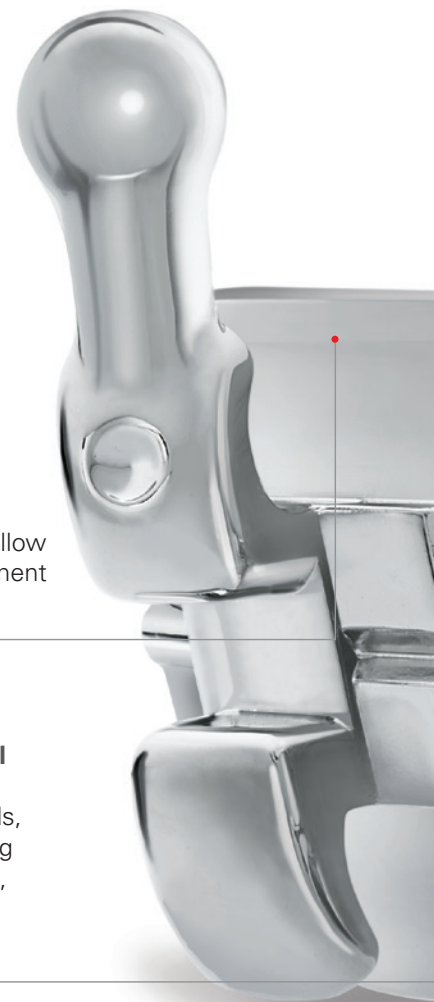

THE *ARIA* BRACKET OFFERS **PATENTED VISUAL CUES** TO ENHANCE PRECISE BRACKET PLACEMENT

Patented Horizontal Scribe Line (HSL) across both the pad surface and the saddle improves alignment, as well as improves the accuracy of bracket height placement.

Rounded mesial and distal entrance minimizes archwire notching thereby reducing "friction"

Smooth and rounded contours allow for maximum patient comfort

Crosshairs for more accurate placement

Bracket height from incisal edge

Correct placement point to line up with the FA Point

KEY POINTS OF THE PATENTED

"HORIZONTAL SCRIBE LINE" (HSL) BY DR. SHIRASUKA



The HSL "Horizontal Scribe Line" of the *Aria* Bracket helps to improve the vertical positioning (occluso-gingival) of the bracket during bonding. In order to ideally position the bracket on the transverse plane (along the long axis of the tooth), we have designed this feature within the bracket base (a perpendicular line from the center of the slot base) to complement bracket height gauges to accurately site on the FA point of the tooth. This reduces improper extrusion or intrusion of teeth, angulation, as well as torque and in-out errors.

The HSL aids in the following:

- Level slot-line up
- Torque direction
- Height positioning with different bracket placement techniques
- Positioning of the pre-molars to help align the marginal ridges
- Removes bracket parallax in higher torque brackets
- Aids in flattening the curve of spee
- Reduces angulation positioning errors

Anatomically contoured base and pad nestles each tooth for positioning that feels right

State-of-the-art pylon base design maximizes surface area, which directly correlates to bond retention

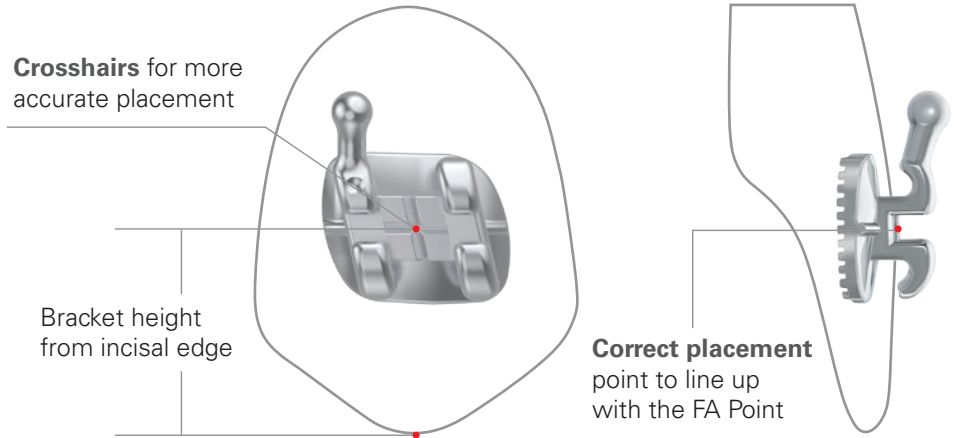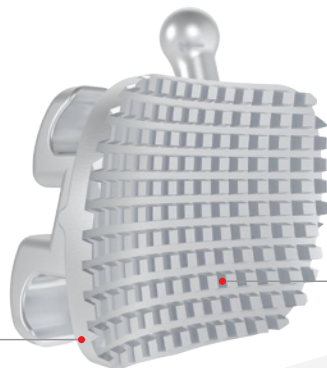

ARIA BRACKETS AND PATIENT SUPPLIES

PRODUCT	ITEM NUMBER
Aria 15X Demo Model	631-716X
Aria Typodont (7x7)	631-716

Aria Brackets Availability:

- .018 slot size, Roth prescription only, and .022 slot size in both Roth and McLaughlin, Bennett, & Trevisi prescriptions
- Hook and non-hook options on cuspids and bicuspid
- All bicuspid offered with or without off-set pads
- Off-set base for upper and lower cuspids
- 17-4 stainless steel construction
- 5x5 case or in packs of 10 brackets



Premium Packaging



Contact your local distributor today
for Low-Profile Brackets.



For more information on our products and educational offerings, please contact us:

OrthoOrganizers.com

In the U.S.: 800.547.2000

Outside the U.S.: +(1) 760 448 8600

Canada: CERUM 800.661.9567

To email an order outside the U.S.:

IntlOrders@OrthoOrganizers.com

